



PLASTIC SLEEVE FOR ASEPTIC SAMPLING SYSTEM 250125 ^(SHOWN)/250225




This installation guide shows the integration of the plastic sleeve to an existing KEOFITT® Aseptic Sampling System.

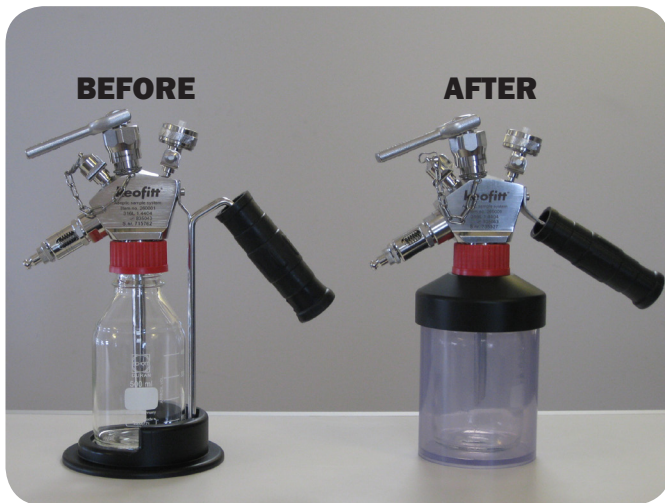

Using the tool supplied unscrew the stainless steel stand and handle as shown.

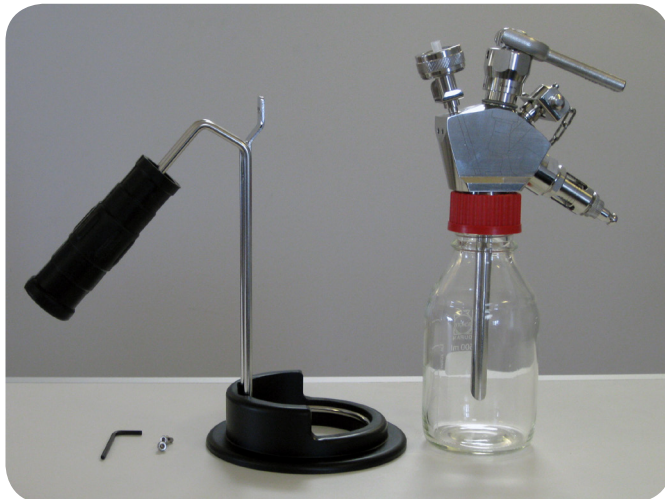

You will not need to reuse the stand, but only the bottle and bottle head shown to your right.

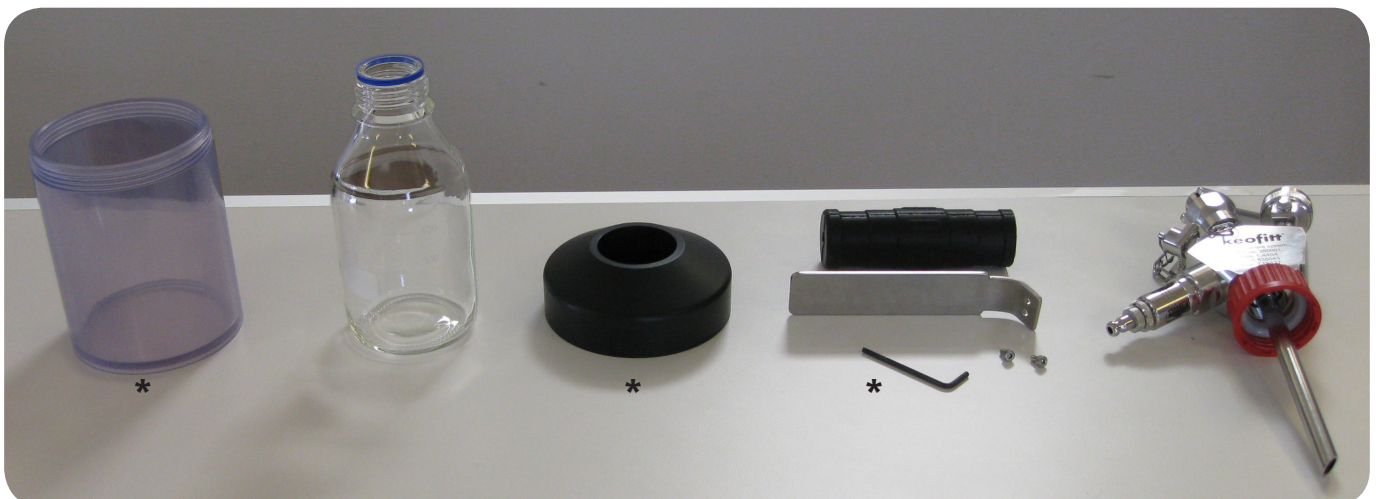

*Item no. 250125 ^(SHOWN)/250225 includes plastic sleeve, plastic top and handle with screws and tool.



Attach the new handle using the screws and the tool supplied.



Mount the plastic top on the bottle neck. The internal thread incorporates separate vertical grooves that allows pressure to exit in case of faulty use.



Mount the bottle with the plastic top onto the screw cap of the Aseptic bottle head as shown. **Please observe that this is the configuration that should be autoclaved. It should not include the plastic sleeve during autoclaving.**



After autoclaving the plastic top and the plastic sleeve should be screwed together.



The Keofitt® Aseptic Sampling System is now ready for use. The plastic sleeve serves two purposes:

- Protection against product and glass.
- Safe and stable stand.

This installation guide only describes the modification of existing Aseptic Sampling Systems.

Please refer to the Aseptic Sampling System User Manual for use of the system.